

### Highlights:

- Results in 5 minutes or less
- Any combination of strips in convenient comb format
- Up to 8 tests - in any combination - from one Common Extraction™:
  - Cry1Ab (Bt)
  - Cry9C (StarLink®)
  - NK 603 (Roundup Ready®)
  - Cry3Bb (YieldGard® Rootworm)
  - Cry1F (Herculex™ I)
  - T25 (LibertyLink®)
  - Cry34 (Herculex RW)
  - mCry3A (Agrisure RW)

### Contents of Kit:

- 4 to 8 QuickStix Strips per comb, packaged 5 combs per foil bag (3 strips per comb or less are packaged 10 combs per canister)
- Sample cups and disposable transfer pipettes

### Items Not Provided:

- Waring blender, model 31BL91 or equivalent
- Glass jar adapter (Eberbach #E8495)
- Glass Mason jars
- Graduated cylinder
- Tap water
- Protective cover for blender jar while grinding



Catalog Number AS 036 TC

## Intended Use

This EnviroLogix QuickComb Kit for bulk grain is designed to extract and detect the presence of certain proteins at the levels typically expressed in genetically modified corn bulk grain. The QuickComb may contain any combination of three to eight of the following QuickStix™:

Protein/Trade Name	Sensitivity
Cry1Ab/Bt11, YieldGard Corn Borer	8 kernels in 800
Cry9C/StarLink	1 kernel in 800
NK 603/Roundup Ready	4 kernels in 800
Cry3Bb/YieldGard Rootworm	4 kernels in 800
Cry1F/Herculex I	4 kernels in 800
T25/LibertyLink	8 kernels in 800
Cry34/Herculex RW	4 kernels in 800
mCry3A/Agrisure RW	8 kernels in 800

## How the Test Works

In order to detect the proteins expressed by genetically modified bulk grain, the sample must first be extracted to solubilize the protein. All QuickStix included in the QuickComb are specially formulated to detect their analytes in a Common Extraction.

Each QuickStix Strip has an absorbent pad at each end. The protective tape with the arrow indicates the end of the strips to insert into the reaction vial. The sample will travel up the membrane strips and be absorbed into the larger pad at the top of the strip. The portion of the strip between the protective tape and the absorbent pad at the top of the strip is used to view the reactions as described under “Interpreting the Results”. Please avoid bending the comb.

## Sample Preparation

### Step 1: Determine Number and Size of Sub-samples

1. Collect a composite sample according to USDA/GIPSA instructions found in the reference documents listed in the margin on Page 2.
2. The following is another helpful reference for use in designing a sampling plan: Remund, K.M., Dixon, D.A., Wright D.L., Holden, L.R. “Statistical considerations in seed purity testing for transgenic traits”, Seed Science Research, June 2001, Vol. 11 No.2, pp. 101-119.
3. To select the appropriate sample size, determine the purity standard and the degree of confidence required. Confidence level means the statistical probability that the % of GMO corn in the lot is below the selected purity standard. This calculation should be done for each trait tested, then choose the largest sub-sample volume result.

*For sampling scenarios for different screening or confidence levels, refer to the USDA/GIPSA Excel spreadsheet described under Step 1 above, or call EnviroLogix Technical Support for assistance.*

*USDA References:*

- <http://archive.gipsa.usda.gov/reference-library/handbooks/grain-insp/grbook1/bk1.pdf> - USDA Grain Inspection Handbook, Book 1, Grain Sampling.
- <http://archive.gipsa.usda.gov/biotech/sample2.htm> - Guidance document entitled Sampling for the Detection of Biotech Grains.
- <http://archive.gipsa.usda.gov/biotech/sample1.htm> - Practical Application of Sampling for the Detection of Biotech Grains.
- <http://archive.gipsa.usda.gov/biotech/samplingplan1.xls> - This website provides a simple to use Sample Planner (29K Excel Spreadsheet). The planner allows you to enter different assumptions in terms of sample size, number of samples, acceptable quality level and to determine the probability of accepting lots with given concentration levels. It also plots the probabilities in graph form for easy interpretation. Specific data can be saved for documentation and future analyses.

**Corn**

*Grams of Corn x 1.5 = mL of water*

*For example:*

*(100 x 0.25) = 25g x 1.5 = 38mL water*

## Step 2: Determine Sub-sample Weight, Jar Size, Grind Times and Water Volume Needed for Sample Preparation

1. Determine the **average weight** of the grain from the lot to be tested. Count and weigh 100 kernels/seeds, then divide by 100.
2. Calculate the sub-sample weight (g) needed for testing, (number of seeds X **average seed weight**).
3. Choose an appropriate jar size for your sample based upon Table 2.
4. Calculate water volume needed for sample preparation. The water volume to sample weight is a ratio of **1.5 to 1**.

*Example Calculation using a 100 kernel sub-sample with an average kernel weight of 0.25g.*  
 $0.25g \times 100 = 25g \times 1.5ml = 38mL$  water for extraction

**Table 2**

Commodity	Sub-sample Weight (g)	Jar Size (oz.)	Grind Time (sec.)
Corn	10-25	4	30
	25-65	8	30
	65-250	32	45

## Step 3: Prepare the Sample


1. Weigh sample into the appropriate size glass Mason jar and attach jar adapter with blade.
2. Place unit on the Waring blender (or equivalent) and cover with protective cover.
3. Grind sample on high speed for specified grinding time or until all whole kernels are broken.
4. Add the volume of tap water calculated by the formula at left.
5. Cap and shake jar vigorously until the entire sample is wet (20-30 seconds, depending on the number of grains). Sample will begin to settle immediately and liquid can be drawn off at that time.
6. Transfer 6 to 12 ml of the liquid portion from above the settled sample into the sample cup. The level should be above the bottom of the arrows but below the top of the lower colored portion. Avoid pulling up particles.
7. To prevent cross-contamination, thoroughly clean blender parts and jars to remove dust and residue prior to preparation of each sample, and use a new sample cup for each. If pipetting, use a new tip or disposable pipette for each sample.

## How to Run the QuickComb Test


1. Remove a QuickComb from the container (foil bag or canister) and return unused combs to original container (avoid handling loose comb end). Use the blank space on the back of the comb to label sample, if desired. Place the comb of strips into the sample cup, being sure to insert the end indicated by the arrows on the protective tape.
2. After inserting the comb into the extract, liquid will travel up the membrane strips toward the absorbent pads at the top of the strips. Soon after complete wetting of the membranes, lines will appear on the membranes approximately 1/4 inch below the top absorbent pad. This is the Control Line.



Transfer extract to cup,  
about 2 pipettefuls

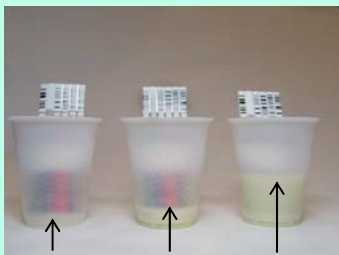

too little extract    just right    too much extract

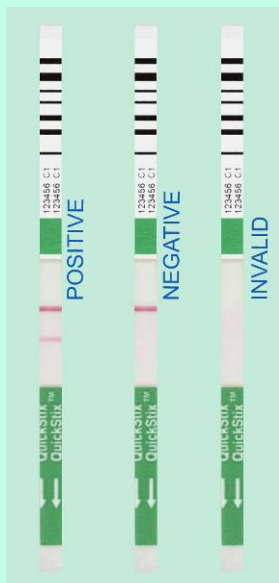


Any clearly discernable pink  
Test Line is considered positive

- The results should develop within 5 minutes. Allow the strips to develop for a full 5 minutes before making final negative assay interpretations. Strongly positive samples may show results much sooner. Remove the QuickComb from the cup to read. To retain the combs, cut off and discard the bottom section of each strip covered by the arrow tape.

## Interpreting the Results

Development of the Control Lines within 5 minutes indicates that the strips have functioned properly. Any strip that does not develop a Control Line should be discarded and the sample re-tested using another strip.

If the extract is from a sample containing at least the detection level of the strip's analyte on the QuickComb, a second line (Test Line) will develop on the membrane strip between the Control Line and the protective tape. *The results should be interpreted as positive for that strip's protein expression.*

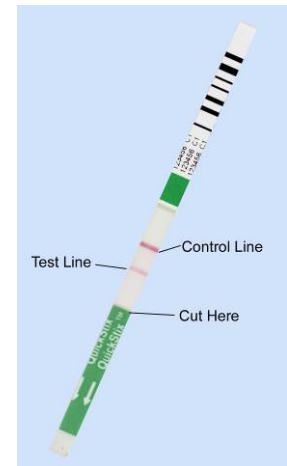

If the extract is from a sample containing less than the listed detection levels, the strip will only develop a Control Line.

## Kit Storage

This QuickComb Kit should be stored at room temperature, or refrigerated for longer shelf life. Please note the shelf life on the kit box for each storage temperature. The kit may be used in field applications; however, prolonged exposure to high temperatures may adversely affect the test results. Important: do not open the container (foil bag or canister) until you are ready to use the combs. Allow container to come to room temperature before opening to prevent condensation. Immediately re-seal unused QuickCombs in the container.

## Precautions and Notes

- This kit is designed to screen for presence or absence only, and is not meant to be quantitative.
- As with all tests, it is recommended that results be confirmed by an alternate method when necessary.
- The assay has been optimized to be used with the protocol provided in the kit. Deviation from this protocol may invalidate the results of the test.
- The results generated through the proper use of this kit reflect the condition of the working sample directly tested. Extrapolation as to the condition of the originating lot, from which the working sample was derived, should be based on sound sampling procedures and statistical calculations which address random sampling effects, non-random seed lot sampling effects and assay system uncertainty. A negative result obtained when properly testing the working sample does not necessarily mean the originating lot is entirely negative for the analyte or protein in question.
- A strong positive result may safely be interpreted in as little as 2 minutes after sample addition. It is not safe, however, to interpret negative results prior to 5 minutes.
- Protect all components from hot or cold extremes of temperature when not in use. Do not leave in direct sunlight or in a vehicle.



**For Technical Support  
Contact Us At:**

**EnviroLogix**  
500 Riverside Industrial  
Parkway  
Portland, ME 04103-1486  
USA

**Tel: (207) 797-0300**  
**Toll Free: 866-408-4597**  
**Fax: (207) 797-7533**

e-mail:  
**[info@envirologix.com](mailto:info@envirologix.com)**

website:  
**[www.envirologix.com](http://www.envirologix.com)**



## LIMITED WARRANTY

EnviroLogix Inc. (“EnviroLogix”) warrants the products sold hereunder (“the Products”) against defects in materials and workmanship when used in accordance with the applicable instructions for a period not to extend beyond a product’s printed expiration date. If the Products do not conform to this Limited Warranty and the customer notifies EnviroLogix in writing of such defects during the warranty period, including an offer by the customer to return the Products to EnviroLogix for evaluation, EnviroLogix will repair or replace, at its option, any product or part thereof that proves defective in materials or workmanship within the warranty period.

**NEITHER ENVIROLOGIX NOR MONSANTO MAKES ANY OTHER WARRANTIES, EXPRESSED OR IMPLIED, INCLUDING BUT NOT LIMITED TO ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE.** The warranty provided herein and the data, specifications and descriptions of EnviroLogix products appearing in EnviroLogix published catalogues and product literatures are EnviroLogix’ sole representations concerning the Products and warranty. No other statements or representations, written or oral, by EnviroLogix’ employees, agents or representatives, except written statements signed by a duly authorized officer of EnviroLogix Inc., are authorized; they should not be relied upon by the customer and are not a part of the contract of sale or of this warranty.

EnviroLogix does not warrant against damages or defects arising in shipping or handling, or out of accident or improper or abnormal use of the Products; against defects in products or components not manufactured by EnviroLogix, or against damages resulting from such non-EnviroLogix made products or components. EnviroLogix passes on to customer the warranty it received (if any) from the maker thereof of such non-EnviroLogix made products or components. This warranty also does not apply to Products to which changes or modifications have been made or attempted by persons other than pursuant to written authorization by EnviroLogix.

**THIS WARRANTY IS EXCLUSIVE.** The sole and exclusive obligation of EnviroLogix shall be to repair or replace the defective Products in the manner and for the period provided above. EnviroLogix shall not have any other obligation with respect to the Products or any part thereof, whether based on contract, tort, strict liability or otherwise. Under no circumstances, whether based on this Limited Warranty or otherwise, shall EnviroLogix be liable for incidental, special, or consequential damages.

This Limited Warranty states the entire obligation of EnviroLogix with respect to the Products. If any part of this Limited Warranty is determined to be void or illegal, the remainder shall remain in full force and effect.

## License

EnviroLogix has developed some tests in this kit using proprietary reagents as well as reagents licensed from Monsanto Company.

*Roundup Ready and YieldGard are registered trademarks of Monsanto Technology, LLC*

*Agrisure is a registered trademark of Syngenta Group Company*

*LibertyLink is a trademark of Bayer*

*Herculex is a trademark of Dow AgroSciences LLC*

*StarLink is a trademark of Aventis CropScience*

*EnviroLogix, the EnviroLogix logo, QuickComb, QuickStix, and Common Extraction are trademarks of EnviroLogix Inc.*

© EnviroLogix 2008