

Highlights:

- Results in 5 minutes or less
- Available as 100-strip kits, or in bulk packaging

Contents of Kit:

- 100 QuickStix Strips packed in two moisture-resistant canisters
- 100 fixed volume transfer pipettes (250 microliters)
- 100 reaction vials

Items Not Provided:

- Oster® Sunbeam blender, 4000 series or equivalent*
- Ice blade with rubber gasket*
- ½ pint (8 oz.) Mini blender jars*
- Threaded bottom cap*
- Protective cover for blender jar while grinding
- Tap water

* Available through EnviroLogix, see "Notes"



Catalog Number AS 017 BG

Intended Use

The EnviroLogix QuickStix Kit for Roundup Ready Canola Bulk Grain is designed to extract and detect the presence of the Roundup Ready protein at the levels typically expressed in genetically modified canola. The sensitivity of this QuickStix Kit is 0.1% (i.e. one Roundup Ready seed in 1000 conventional canola seeds). For detection of Roundup Ready in canola plant tissues and individual seeds, please use QuickStix Cat# AS 017 LS.

How the Test Works

In order to detect the Roundup Ready protein with this QuickStix Kit, the sample must first be crushed and extracted with tap water to solubilize the protein.

Each QuickStix Strip has an absorbent pad at each end. The protective tape with the arrow indicates the end of the strip to insert into the reaction vial. The sample will travel up the membrane strip and be absorbed into the larger pad at the top of the strip. The portion of the strip between the protective tape and the absorbent pad at the top of the strip is used to view the reactions as described under "Interpreting the Results".

Sample Preparation

Step 1: Determine Number and Size of Sub-samples

1. Collect a composite sample according to USDA/ GIPSA instructions found in the following reference documents:
 - <http://archive.gipsa.usda.gov/reference-library/handbooks/grain-insp/grbook1/bk1.pdf> - USDA Grain Inspection Handbook, Book 1, Grain Sampling.
 - <http://archive.gipsa.usda.gov/biotech/sample2.htm> - Guidance document: "Sampling for the Detection of Biotech Grains".
 - <http://archive.gipsa.usda.gov/biotech/sample1.htm> - Practical Application of Sampling for the Detection of Biotech Grains.
 - <http://archive.gipsa.usda.gov/biotech/samplingplan1.xls> - This website provides a simple to use Sample Planner (29K Excel Spreadsheet). The planner allows entry of different assumptions in terms of sample size, number of samples, acceptable quality level and to determine the probability of accepting lots with given concentration levels. It also plots the probabilities in graph form for easy interpretation. Specific data can be saved for documentation and future analyses.
2. The following is a helpful reference for use in designing a sampling plan: Remund, K.M., Dixon, D.A., Wright D.L., Holden, L.R. "Statistical considerations in seed purity testing for transgenic traits", Seed Science Research, June 2001, Vol. 11 No.2, pp. 101-119.
3. To select the appropriate sample size, first determine the purity standard and the degree of confidence required. Confidence level means the statistical probability that the % of Roundup Ready in the lot is below the selected purity standard. Table 1 provides a guideline for determining the number of sub-samples necessary



To calculate the average seed weight:

- Count & weigh 100 seeds
- Divide total weight by 100



Weigh sample into blender jar, grind



Add calculated volume of water

Buffer Volume Formula

$$\text{Grams canola seed} \times 5 = \text{mL of tap water to add}$$

For example: a 332-seed sample with an average weight of 0.0033g:
 $(332 \times 0.0033 = 1.06 \text{ g} \times 5 = 5.3 \text{ mL tap water})$



Shake to mix, allow to settle briefly

to provide effective screening for different GM concentrations at the 95% and 99% confidence levels.

Table 1 – Canola - Number of seeds required

Testing 1 sub-sample with that sub-sample being negative (-)

Confidence Level (%)	Roundup Ready Screening Level		
	5%	3%	0.9%
95%	59	99	332
99%	90	152	510

Testing 2 sub-samples with 1 or both being negative (-/-) or (-/+)

	5%	3%	0.9%
95%	72	120	407
99%	104	174	586

Testing 3 sub-samples with 2 or 3 sub-samples being negative (-/-/-) or (-/-/+)

	5%	3%	0.9%
95%	39	66	221
99%	55	93	315

Note: Screening at the 0.1% Roundup Ready concentration level, with 95% confidence, requires 3 sub-samples of 1000 seeds with all sub-samples being negative (-/-/-).

For other sampling scenarios or different screening or confidence levels, refer to the USDA/GIPSA spreadsheet described under Step 1 above, or call EnviroLogix Technical Support for assistance.

Step 2: Determine Sub-sample Weight

1. Determine average weight of individual grain to be tested (see formula, above left).
2. Calculate the sub-sample weight (g) needed for testing, (number of seeds X average seed weight).

Step 3: Prepare the Sample


1. Weigh sample into the blender jar.
2. Put rubber gasket then blade atop the jar.
3. Attach black blender base to jar and hand tighten.
4. Invert and put the blender jar on the blender motor, turn blender on high for 30 seconds.
5. Remove the blender jar, invert and tap. Remove the base and blade assembly. Visually inspect to ensure that every seed has been broken. Add the volume of tap water calculated at left. For sample sizes of fewer than 200 seeds, use a minimum volume of 3.0 mL water to extract.
6. Put white lid on the blender jar.
7. Shake the sample 30 times.



Pipette sample into vial



Add QuickStix and read the results

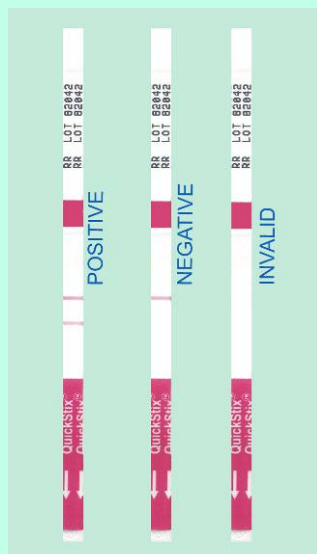

Any clearly discernable pink Test Line is considered positive

- Let sample settle (approximately 30 seconds). Dispense 250 μ L extract into reaction vial using disposable fixed volume pipette provided with the kit; avoid particulates.
- To prevent cross-contamination, use a new transfer pipette and reaction vial for each sample.

How to Run the QuickStix Strip Test

- Allow refrigerated canisters to come to room temperature before opening. Remove the QuickStix Strips to be used. Avoid bending the strips. Reseal the canister immediately.
- Place the strip into the reaction vial. The sample will travel up the strip. Reaction vials will stand on their own or may be inserted into the cardboard racks provided.
- Allow the strip to develop for 5 minutes before making final assay interpretations. Positive sample results may become obvious much more quickly.
- To retain the strip, cut off and discard the bottom section of the strip covered by the arrow tape.


NOTE: Use extreme caution to prevent sample-to-sample cross-contamination with grain, fluids, or disposables.

Interpreting the Results

Development of the Control Line within 5 minutes indicates that the strip has functioned properly. Any strip that does not develop a Control Line should be discarded, and the sample re-tested using another strip.

If the extract is from a sample containing at least 0.1% Roundup Ready-modified canola, a second line (Test Line) will develop on the membrane strip between the Control Line and the protective arrow tape. *The results should be interpreted as positive for the presence of Roundup Ready protein.*

If the extract is from a negative sample, the strip will only show the control line.

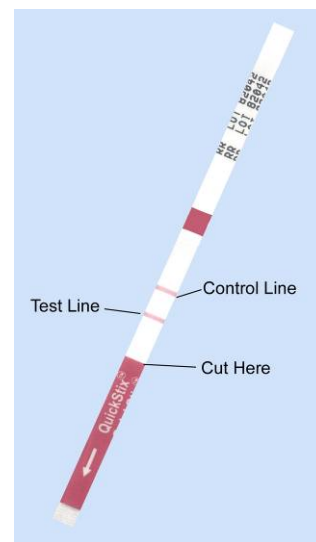

Kit Storage

QuickStix can be stored at room temperature, or refrigerated for a longer shelf life. Note the shelf life on the kit box for each storage temperature. The kit may be used in field applications; however, prolonged exposure to high temperatures may adversely affect the test results. Do not open the desiccated canister until ready to use the test strips.

Precautions and Notes

- This kit is designed to screen for presence or absence only, and is not meant to be quantitative.
- As with all tests, it is recommended that results be confirmed results by an alternate method if necessary.
- The assay has been optimized to be used with the protocol provided in the kit. Deviation from this protocol may invalidate the results of the test.
- The results generated through the proper use of this diagnostic tool reflect the condition of the working sample directly tested. Extrapolation as to the condition



of the originating lot, from which the working sample was derived, should be based on sound sampling procedures and statistical calculations which address random sampling effects, non-random seed lot sampling effects and assay system uncertainty. A negative result obtained when properly testing the working sample does not necessarily mean the originating lot is entirely negative for the analyte or protein in question.

- A strong positive result may safely be interpreted in as little as 5 minutes after sample addition. It is not safe to interpret weak positive or negative results prior to 10 minutes.
- DO NOT leave in direct sunlight or in vehicle. Protect all components from hot or cold extremes of temperature when not in use.

Accessory Items Available

- Oster Sunbeam Blender Model 4094 (Cat. # ACC 044)
- 8 oz. plastic blender jars with caps for Oster Sunbeam (10/pk, Cat. # ACC 046)
- Replacement blade for Oster blender (Cat. # ACC 045)





**For Technical Support
Contact Us At:**

EnviroLogix
500 Riverside Industrial
Parkway
Portland, ME 04103-
1486 USA
Tel: (207) 797-0300
Toll Free: 866-408-4597
Fax: (207) 797-7533

e-mail:
info@envirologix.com

website:
www.envirologix.com



LIMITED WARRANTY

EnviroLogix Inc. ("EnviroLogix") warrants the products sold hereunder ("the Products") against defects in materials and workmanship when used in accordance with the applicable instructions for a period not to extend beyond a product's printed expiration date. If the Products do not conform to this Limited Warranty and the customer notifies EnviroLogix in writing of such defects during the warranty period, including an offer by the customer to return the Products to EnviroLogix for evaluation, EnviroLogix will repair or replace, at its option, any product or part thereof that proves defective in materials or workmanship within the warranty period.

NEITHER ENVIROLOGIX NOR MONSANTO MAKES ANY OTHER WARRANTIES, EXPRESSED OR IMPLIED, INCLUDING BUT NOT LIMITED TO ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. The warranty provided herein and the data, specifications and descriptions of EnviroLogix products appearing in EnviroLogix published catalogues and product literature are EnviroLogix' sole representations concerning the Products and warranty. No other statements or representations, written or oral, by EnviroLogix' employees, agents or representatives, except written statements signed by a duly authorized officer of EnviroLogix Inc., are authorized; they should not be relied upon by the customer and are not a part of the contract of sale or of this warranty.

EnviroLogix does not warrant against damages or defects arising in shipping or handling, or out of accident or improper or abnormal use of the Products; against defects in products or components not manufactured by EnviroLogix, or against damages resulting from such non-EnviroLogix made products or components. EnviroLogix passes on to customer the warranty it received (if any) from the maker thereof of such non-EnviroLogix made products or components. This warranty also does not apply to Products to which changes or modifications have been made or attempted by persons other than pursuant to written authorization by EnviroLogix.

THIS WARRANTY IS EXCLUSIVE. The sole and exclusive obligation of EnviroLogix shall be to repair or replace the defective Products in the manner and for the period provided above. EnviroLogix shall not have any other obligation with respect to the Products or any part thereof, whether based on contract, tort, strict liability or otherwise. Under no circumstances, whether based on this Limited Warranty or otherwise, shall EnviroLogix be liable for incidental, special, or consequential damages.

This Limited Warranty states the entire obligation of EnviroLogix with respect to the Products. If any part of this Limited Warranty is determined to be void or illegal, the remainder shall remain in full force and effect.

License

EnviroLogix has developed this kit using proprietary reagents as well as reagents licensed from Monsanto Company.

*Roundup Ready is a registered trademark of Monsanto Technology, LLC
EnviroLogix, the EnviroLogix logo, and QuickStix are trademarks of EnviroLogix Inc.*

© Envirologix 2009