

### Contents of Kit:

- 100 QuickStix Strips packed in two moisture-resistant canisters
- Sample cups and disposable transfer pipettes

### Items Not Provided:

- Bunn grinder or equivalent
- Graduated cylinder
- Tap water
- QuickScan System

Catalog Number AQ 074 BG

## Intended Use

The QuickStix Kit for QuickScan – eCry3.1Ab detects and quantifies the eCry3.1Ab protein at the levels typically expressed in Event 5307 (Agrisure Duracade™) corn grain. Corn derived from transformation Event 5307 contains the gene *ecry3.1Ab* which encodes eCry3.1Ab protein expression in the plant tissues. This protein is an engineered chimera of modified Cry3A (mCry3A) and Cry1Ab proteins. The sensitivity of these QuickStix Strips is 0.1% Event 5307 corn (i.e. one kernel in 1000).

## How the Test Works


In order to detect the eCry3.1Ab protein expressed by Duracade corn, the sample must first be extracted to solubilize the protein. Each QuickStix Strip has an absorbent pad at each end. The protective tape with the arrow indicates the end of the strip to insert into the sample cup. The sample will travel up the membrane strip and be absorbed into the larger pad at the top of the strip. The portion of the strip between the protective tape and the absorbent pad at the top of the strip is used to interpret the reactions as described under “Interpreting the Results.” Results are scanned and interpreted quantitatively with the EnviroLogix QuickScan System. Please avoid bending the strips.

## Sample Preparation

1. Collect a composite sample according to USDA/GIPSA instructions found in the reference documents listed in Precautions and Notes.
2. Determine the **average weight** of the grain from the lot to be tested. Count and weigh 100 kernels/seeds, then divide by 100. Corn kernels vary widely from region to region and from year to year. Be sure to test the desired number of kernels by checking average kernel weight periodically.
3. Calculate the sub-sample weight (g) needed for testing, (number of seeds X **average seed weight**). A sub-sample size of at least 200 grams is required to create enough extract for the test. Note: it is necessary to grind more than the calculated sub-sample weight—subsample is weighed out after grinding (Step 5).
4. Calculate water volume needed for sample preparation. The Common Extraction Method calls for a water volume to sample weight ratio of **1.5 to 1** (see example, right).
5. Grind corn using the Auto-Drip setting on the Bunn grinder (or equivalent), then weigh out the subsample calculated in Step 3. The sample should be the consistency of coffee grounds – 60-70% of the sample should pass through a 20-mesh sieve.
6. Weigh out subsample calculated in Step 3, then place sub-sample into an appropriately sized jar or stand-up, zip-type plastic bag and add the volume of tap water calculated in Step 4.
7. Close the jar or bag and shake vigorously for 30 seconds, then allow sample to settle for another 30 seconds.
8. Assemble cardboard holder and lift the support. (Note: save and reuse holders.) On a flat surface, insert sample cup in the space provided, then transfer approximately 35 mL of the liquid portion from above the settled sample into the sample cup. Use a fresh pipette from the kit, or pour slowly and carefully, until the depth of the sample is level with the upper surface of the cardboard holder. Important: Avoid transferring particles as much as possible, and after transfer, allow the liquid in the sample cup to settle for 30 seconds so that any particles will settle at the bottom of the cup.
9. To prevent cross-contamination, thoroughly clean blender parts and jars to remove dust and residue prior to preparation of each sample, and use a new sample cup for each. If pipetting, always use a new pipette for each sample.

Corn Common Extraction™  
g of corn X 1.5 = mL of water

*Example:*  
800 kernel sub-sample with an  
average kernel weight of 0.3g.  
 $0.3 \text{ g} \times 800 = 240 \text{ g} \dots$   
 $240 \times 1.5 = 360 \text{ mL water}$

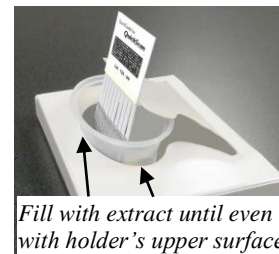

## How to Run the QuickStix Strip Test

1. Allow refrigerated canisters to come to room temperature before opening. Remove the QuickStix Strips to be used. Avoid bending the strips. Reseal the canister immediately.
2. Place the strip into the sample cup provided, being sure to insert the end indicated by the arrows on the protective tape. The sample will travel up the strip.
3. Allow the strip to develop for 5 minutes before making final assay interpretations.

4. Immediately cut off and discard the bottom section of the strip covered by the arrow tape and place in QuickScan Reader. Strips must be read while still wet.

**NOTE: Use extreme caution to prevent sample-to-sample cross-contamination with grain, fluids, or disposables.**

## Interpreting the Results

Development of the Control Lines within 5 minutes indicates that the strips have functioned properly. Any strip that does not develop a Control Line should be discarded and the sample re-tested using another device.

Results are scanned and interpreted quantitatively with the QuickScan System. Place QuickStix Strip into the carrier, slide in, and press “Read Test” on the screen. QuickScan will return a result as “% GMO” or “<LOD” (less than the Limit of Detection). Please consult the QuickScan User Manual for details.

## Kit Storage

QuickStix can be stored at room temperature, or refrigerated for a longer shelf life. Note the shelf life on the kit box for each storage temperature. The kit may be used in field applications; however, prolonged exposure to high temperatures may adversely affect the test results. Do not open the desiccated canister until ready to use the test strips.

## Precautions and Notes

- This kit is designed to give a quantitative result using the QuickScan System and is not intended to be visually interpreted.
- As with all tests, it is recommended that results be confirmed by an alternate method when necessary.
- The assay has been optimized to be used with the protocol provided in the kit. Deviation from this protocol may invalidate the results of the test.
- The results generated through the proper use of this kit reflect the condition of the working sample directly tested. Extrapolation as to the condition of the originating lot, from which the working sample was derived, should be based on sound sampling procedures and statistical calculations which address random sampling effects, non-random seed lot sampling effects and assay system uncertainty. A negative result obtained when properly testing the working sample does not necessarily mean the originating lot is entirely negative for the analyte or protein in question.
- Protect all components from hot or cold extremes of temperature when not in use. Do not leave in direct sunlight or in a vehicle.
- USDA References:
  - [www.gipsa.usda.gov/fgis/handbook/BK1/Book1\\_2015-09-18.pdf](http://www.gipsa.usda.gov/fgis/handbook/BK1/Book1_2015-09-18.pdf) - USDA Grain Inspection Handbook, Book 1, Grain Sampling.
  - [www.gipsa.usda.gov/fgis/biotech/sample2.htm](http://www.gipsa.usda.gov/fgis/biotech/sample2.htm) - Guidance document entitled Sampling for the Detection of Biotech Grains.
  - [www.gipsa.usda.gov/fgis/biotech/sample1.htm](http://www.gipsa.usda.gov/fgis/biotech/sample1.htm) - Practical Application of Sampling for the Detection of Biotech Grains.
  - [www.gipsa.usda.gov/fgis/biotech/samplingplan1.xlsx](http://www.gipsa.usda.gov/fgis/biotech/samplingplan1.xlsx) - This website provides a simple to use Sample Planner (29K Excel Spreadsheet). The planner allows you to enter different assumptions in terms of sample size, number of samples, acceptable quality level and to determine the probability of accepting lots with given concentration levels. It also plots the probabilities in graph form for easy interpretation. Specific data can be saved for documentation and future analyses.

**LIMITED WARRANTY**

EnviroLogix Inc. (“EnviroLogix”) warrants the products sold hereunder (“the Products”) against defects in materials and workmanship when used in accordance with the applicable instructions for a period not to extend beyond a product’s printed expiration date. If the Products do not conform to this Limited Warranty and the customer notifies EnviroLogix in writing of such defects during the warranty period, including an offer by the customer to return the Products to EnviroLogix for evaluation, EnviroLogix will repair or replace, at its option, any product or part thereof that proves defective in materials or workmanship within the warranty period.

**ENVIROLOGIX MAKES NO OTHER WARRANTIES, EXPRESS OR IMPLIED, INCLUDING BUT NOT LIMITED TO ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE.** The warranty provided herein and the data, specifications and descriptions of EnviroLogix products appearing in EnviroLogix published catalogues and product literature are EnviroLogix’ sole representations concerning the Products and warranty. No other statements or representations, written or oral, by EnviroLogix’ employees, agents or representatives, except written statements signed by a duly authorized officer of EnviroLogix Inc., are authorized; they should not be relied upon by the customer and are not a part of the contract of sale or of this warranty.

EnviroLogix does not warrant against damages or defects arising in shipping or handling, or out of accident or improper or abnormal use of the Products; against defects in products or components not manufactured by EnviroLogix, or against damages resulting from such non-EnviroLogix made products or components. EnviroLogix passes on to customer the warranty it received (if any) from the maker thereof of such non-EnviroLogix made products or components. This warranty also does not apply to Products to which changes or modifications have been made or attempted by persons other than pursuant to written authorization by EnviroLogix.

**THIS WARRANTY IS EXCLUSIVE.** The sole and exclusive obligation of EnviroLogix shall be to repair or replace the defective Products in the manner and for the period provided above. EnviroLogix shall not have any other obligation with respect to the Products or any part thereof, whether based on contract, tort, strict liability or otherwise. Under no circumstances, whether based on this Limited Warranty or otherwise, shall EnviroLogix be liable for incidental, special, or consequential damages.

This Limited Warranty states the entire obligation of EnviroLogix with respect to the Products. If any part of this Limited Warranty is determined to be void or illegal, the remainder shall remain in full force and effect.

For Technical Support Contact Us At:

EnviroLogix  
500 Riverside Industrial Parkway  
Portland, ME 04103-1486 USA

Tel: (207) 797-0300 Toll Free (U.S.): 866-408-4597

Fax: (207) 797-7533

e-mail: [info@envirologix.com](mailto:info@envirologix.com)

website: [www.envirologix.com](http://www.envirologix.com)

*Agrisure Duracade™ is a trademark of a Syngenta Group Company.*

*EnviroLogix, the EnviroLogix logo, QuickComb, QuickStix, Common Extraction, and QuickScan are trademarks of EnviroLogix Inc.*

© EnviroLogix 2019